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## WADE- / FUSSMASSAGEGERÄT

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## FOOT MASSAGE DEVICE

Releass

**Thank you for purchasing the Relness foot massage device from MunichSpring.  
Read these user manuals carefully and keep them.**

### Important safeguards

- This device can be used by children aged 12 years and over and by persons with limited physical, sensory or mental capabilities provided they have been given supervision or instructions on how to use the device safely and understand the resulting dangers.
- Please consult your doctor before using if you are taking medication or have any illnesses. You should not use the product before you consult your doctor if you have any of the following problems: pacemaker, recent surgery, epilepsy, recent joint replacements, varicose veins, open wounds, bruises, phlebitis, metal pins or plates, or any physical health concerns.
- Children are not allowed to play with the device.
- If you feel uncomfortable using this device, stop using it and consult a doctor. Ask your doctor if you have any therapeutic questions about using the device. This is a device for household use and for personal use.
- Do not use the device in wet or dusty areas, as this could result in malfunction or electric shock.
- Do not use the device in the bathroom.
- Do not expose the device to heating or direct sunlight.
- Only use accessories that are supplied/ recommended by MunichSpring.
- Do not operate this device if a cord or plug is damaged. Please send the device back to the service center for testing and repair.
- Do not try to repair the device yourself. There are no parts that the user can repair himself.

### Attention





- Do not fall asleep while the device is on.
- Do not use this product under the influence of alcohol or other drugs or if you are physically unwell.
- Moving the product from a cold storage location to a warmer environment can cause malfunction due to water vapor. In this case, let the device acclimatise for an hour in the warmer environment.
- If you find that any component is damaged, or it is in any unusual condition, do not use this product.
- If the power goes out during your massage, do not reconnect the power to avoid further damage.
- Do not use the product if there is thunder and lightning.

### Guarantee Claims

- The guarantee is not transferable.
- This guarantee only applies if the product was purchased directly from MunichSpring.
- Proof of purchase (original receipt) is required for all warranty repairs.
- Any implied warranties are limited to the legally stipulated warranty period as well as to the suitability and general usability.
- This guarantee does not cover damage due to misuse, improper use or unauthorized changes and repairs.
- This guarantee does not apply to accessories that are not part of the massage unit. This guarantee does not apply to rental or business, commercial or institutional use, which does not fall under private use.
- All service work covered by this guarantee must be approved by MunichSpring and carried out only by authorized technicians.
- If no spare parts are available for defective devices, MunichSpring reserves the right to provide replacements instead of repairs or exchanges.

**Main unit & Accessories**

- A 5 control buttons to different functions
- B Vibration area
- C Front sole area rolling and massage
- D Arch area rolling + massage + heat
- E Heel area rolling and massage
- F Air bags
- G Support, could provide 360 degree adjustable support
- H Support adjust button
- I Magnetic area for remote
- J Removable & washable foot pad

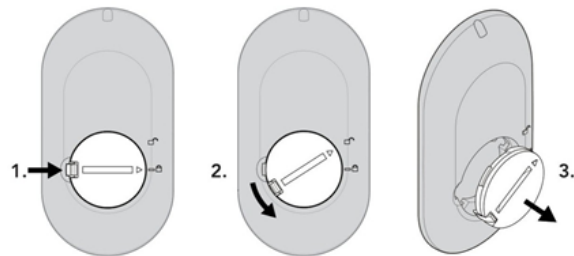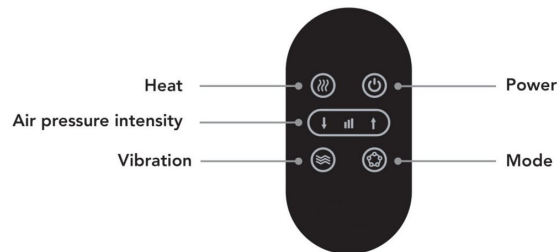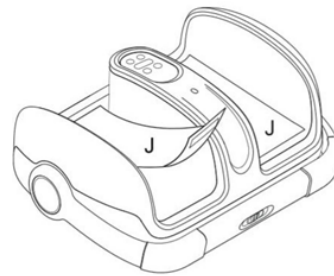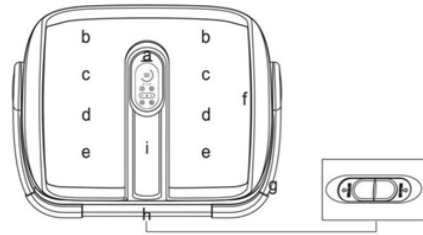Push in clasp, then rotate battery cover to open the battery cover

**Effect and application**

- Power**  
On/Off
- Heat**  
On/Off
- Vibration**  
Off/Low/High
- Air pressure intensity**  
P1 Lowest – P5 Highest
- Mode**
  - Shiatsu + Air Pressure + Vibration
  - Shiatsu only
  - Air pressure only

- Stimulates muscle regeneration in calves, ankles and feet
- Improves blood circulation
- Relieves foot pain caused by plantar fasciitis and heel spurs
- Reduces inflammation
- Relieves tension
- Trigger point and foot reflexology massage
- Alleviates diabetic neuropathy

## Product maintenance

### Storage

- Store the massager in its box or in a safe, dry and cool place. Avoid contact with sharp edges or pointed objects that could damage the surface.
- Please wind up the power cord and stow it under the lower frame of the device. Do not hang the device by the power cord.

### Cleaning

- Do not immerse the device in liquids during cleaning.
- Pull the plug out of the socket and let the device cool down before cleaning. Clean only with a soft, slightly damp cloth.
- Do not let water or other liquids get into the device.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, glass cleaner, furniture polish, or paint thinner to clean the product.
- The fabric cover can be hand washed or washed at 30 °C

## Technical specifications

- Product Name: Relness
- Model Number: RE-214
- Dimensions: 60 x 45,7 x 30,5 cm
- Net weight: : ca. 8,6 kg
- Gross weight: ca. 9,9 kg
- Material: ABS, Metal, Fabric
- Color: gray / black
- Timer: automatic deactivation after approximately 15 min.
- Power: 60W
- Frequency: 50/60 Hz
- Operating voltage: 100 - 240 V

## Environmental Protection

Due to the substances used in electrical and electronic devices that are harmful to the environment and health, devices marked with this symbol must not be disposed of with normal household waste, but must be reprocessed, reused or recycled.

# Beyond relaxation

Your well-being is our priority, and we aim to effectively relieve your pain and to offer you relaxation anytime, anywhere.

With the help of the latest technology and the knowledge of our experts about traditional Asian medicine, we have developed massage and wellness products that accompany you in your stressful everyday life and provide you with comfort and positive energy. Here, we use the power of nature and combine it with our products, to bring your body and mind into balance in order to find your way back to yourself.

With the help of our products, you can relax and relieve your pain and will have more time and energy throughout the day.

We have been on the market since 2013 and for years, with our products, have been able to help countless customers get rid of their pain. We can also be found at many trade fairs, which enables us to actively listen to our customers and collect valuable feedback. As a result, we have not only built trust, but were also able to continuously develop our products by adapting them to the needs and wishes of our customers.

Convince yourself of our products and if you have any suggestions, we always have an open ear for improvement. If you would like to learn more about us and our products, you can contact us at any time or even visit us personally on site and at trade fairs.